

Reflux Sentinel — OLED Menu Reference

Display: SSD1306 / SH1106 1.3" OLED · 128×64 px · Monochrome white-on-black

Controls: Rotary encoder (scroll / click) top-right · Tactile button (back / action) top-left

Mount: Back face of device, visible through PETG clip window

Legend

◀	= back button (top-left)
▶	= encoder click (confirm)
↑	= encoder scroll (navigate)
▶	= selected / active item
▒	= progress / fill bar

1 · First Boot — Onboarding Flow

1.1 Welcome Screen

▒▒▒▒▒▒	▒▒▒▒▒▒	▒▒▒▒▒▒
REFLUX	SENTINEL	
Wearable GERD Monitor		
▶ Press encoder to start		

1.2 Wear-Side Selection

SETUP: Wear Side

► Left side
Right side

Clip faces outward

1.3 Calibration Prompt

SETUP: Calibrate Posture

Stand upright & still

► Click to capture
baseline posture

1.4 Calibration Complete

✓ Calibrated!


Baseline saved

Pitch: 0.0°

Starting in 3s...

1.5 Hotspot QR Code

Wi-Fi Setup

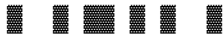


Scan to
connect

Sentinel-AP / 192.168.1.1


QR encodes SSID + WPA2 password for one-tap phone join.

2 · Normal Operation — Watch Face (Home)


2.1 Default Watch Face

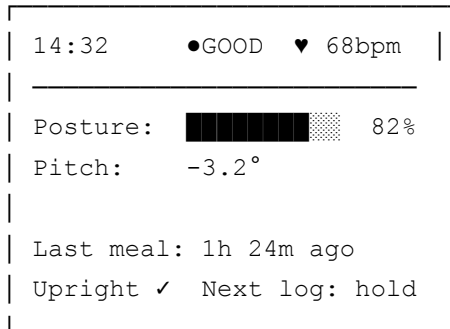

2.2 Watch Face — Alert State

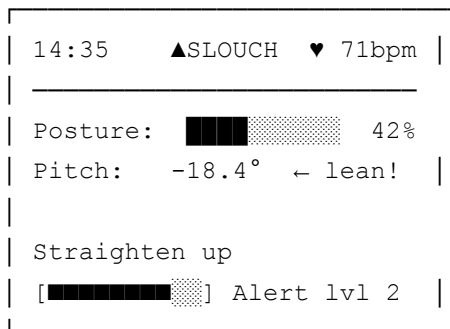


2.3 Watch Face — Post-Meal Countdown

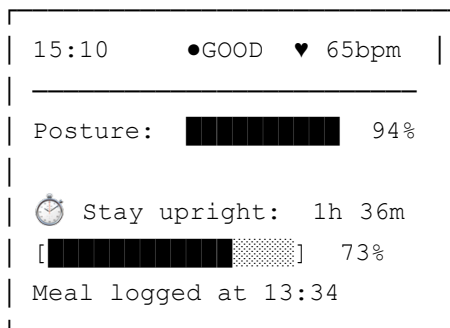


3 • Main Menu

Reached by a single encoder click from the watch face.

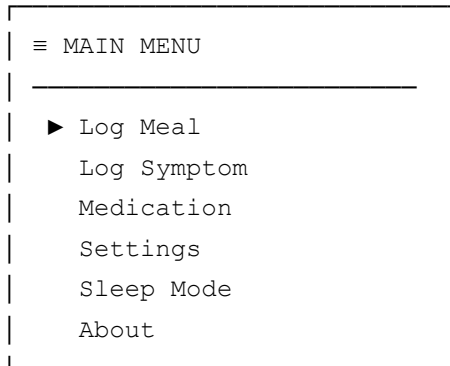


4 • Meal Logging Flow


4.1 Log Meal — Confirm

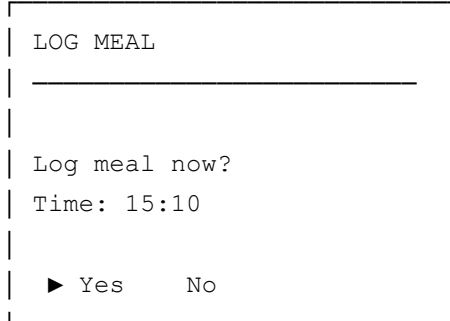


4.2 Food Photo Capture

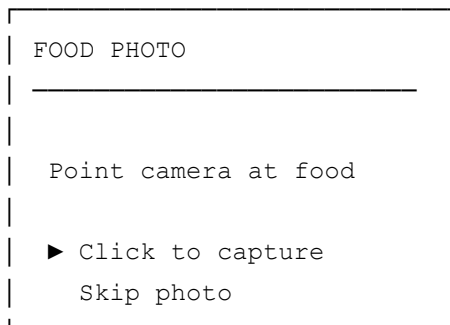


4.3 Analysing Food

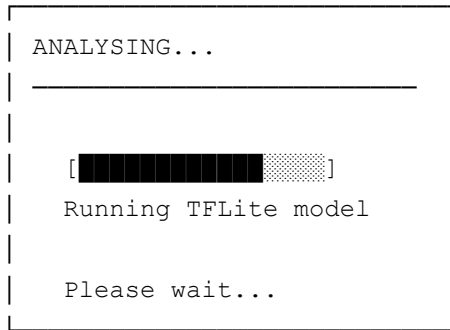


4.4 Food Result — Low Risk

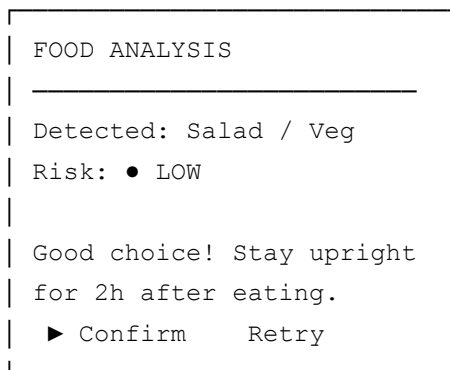

4.5 Food Result — High Risk

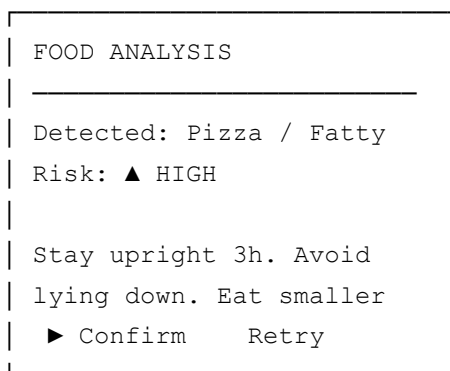

4.6 Meal Saved

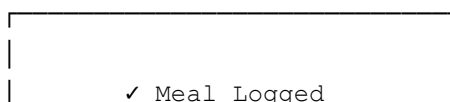

Upright timer started:
Stay upright for 2h 30m

Returning to watch face

5 · Symptom Logging

5.1 Symptom Severity

LOG SYMPTOM

Severity?

► 1 - Mild

2 - Moderate

3 - Severe

5.2 Symptom Type

SYMPTOM TYPE

► Heartburn

Regurgitation

Bloating

Chest pain

Other

5.3 Symptom Saved

✓ Symptom Logged

Severity: 2 - Moderate

Type: Heartburn

Time: 15:14

6 • Medication

6.1 Medication Reminder (auto popup)

 MEDICATION REMINDER

Omeprazole 20mg
Scheduled: 20:00

▶ Tap button to confirm

6.2 Medication Acknowledged

✓ Med Acknowledged

Omeprazole 20mg
Taken at: 20:03

6.3 Medication List (from menu)

MEDICATIONS

▶ Omeprazole 20:00 ✓


Gaviscon 22:00 ...

(Configure on phone)


◀ Back

7 • Settings Menu


7.1 Settings Root

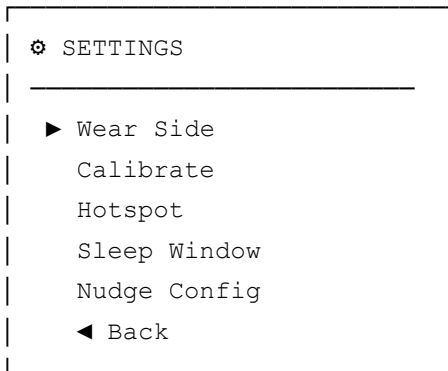


7.2 Wear Side

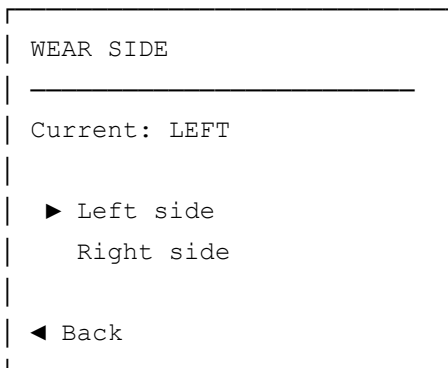


7.3 Recalibrate Posture — Prompt

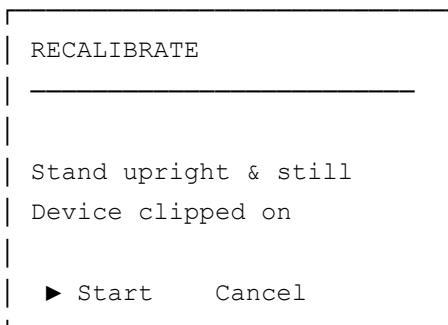


7.3b Recalibrate — Countdown

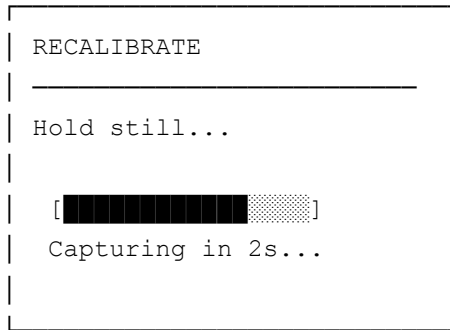


7.3c Recalibrate — Done

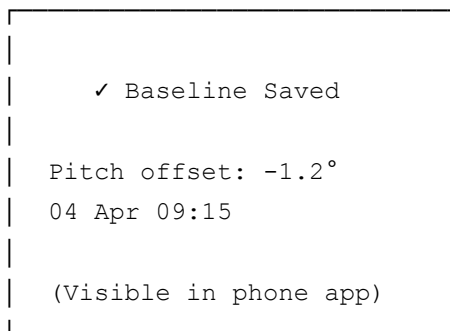

7.4 Hotspot Settings — ON

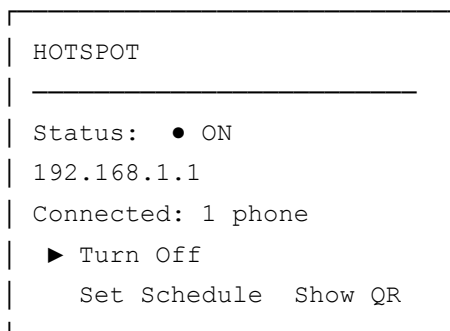

7.4b Hotspot Settings — OFF

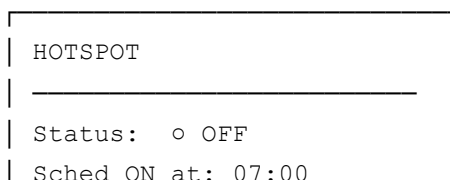

► Turn On Now

Set Schedule

7.4c Hotspot Toggled (transient overlay)

 Hotspot is ON

192.168.1.1

Sentinel-AP

(auto-dismiss 3s)

7.5 Hotspot Schedule

HOTSPOT SCHEDULE

► Mode: Scheduled

ON at: 07:00

OFF at: 23:00

► Click to edit · ◀ Back

7.5b Schedule Edit (encoder scrolls hour/min)

SET ON TIME

[07] : [00]

↑ ↑

scroll scroll

► Confirm ◀ Cancel

7.6 Sleep Window

SLEEP WINDOW
Auto-detect: ON
From: 23:00
Until: 07:00
(Edit on phone PWA)
◀ Back

7.7 Nudge Config


NUDGE CONFIG
Day haptics: • ON
Night haptics: • ON (max 3)
Voice alerts: • ON
(Edit thresholds on phone)
◀ Back

8 • Sleep Mode


8.1 Sleep Mode — Auto Detect

 SLEEP MODE
Detected: lying flat
Roll: +12° (Right side)
Position: ▲ NOT IDEAL
Try: left side sleep


8.2 Sleep Mode — Good Position

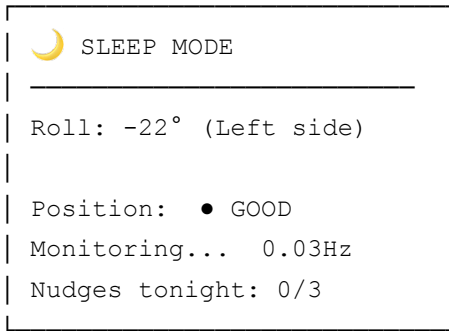

8.3 Sleep Mode — Grace Period

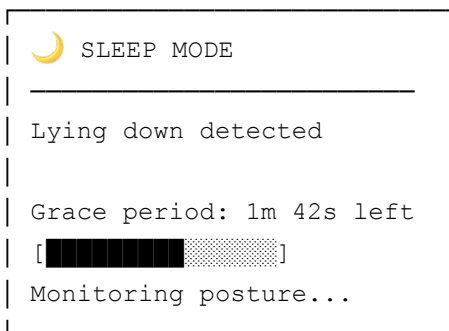


9 · About Screen

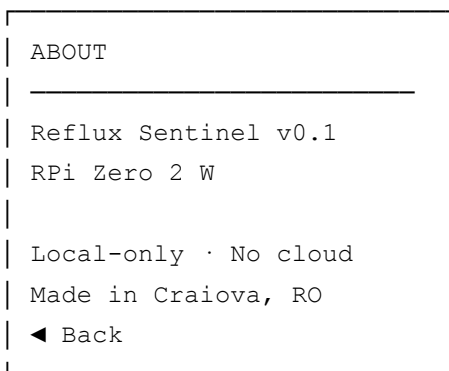


10 · Transient Overlays (non-navigable)


10.1 Hydration Reminder

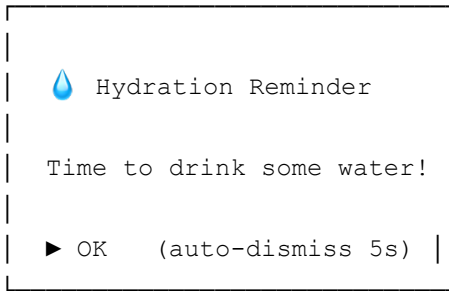


10.2 Hotspot Connected

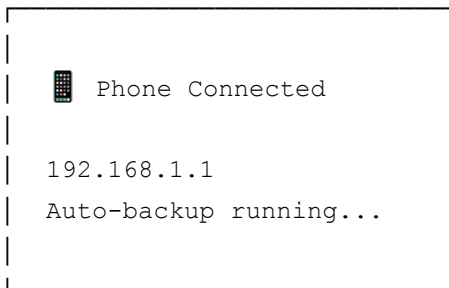


10.3 Low Power Warning

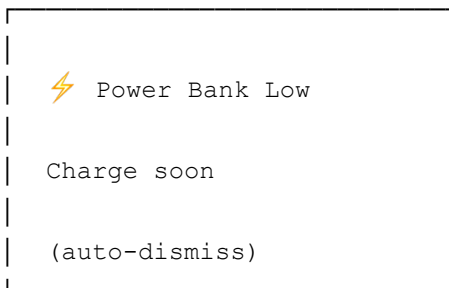

OLED Navigation Map

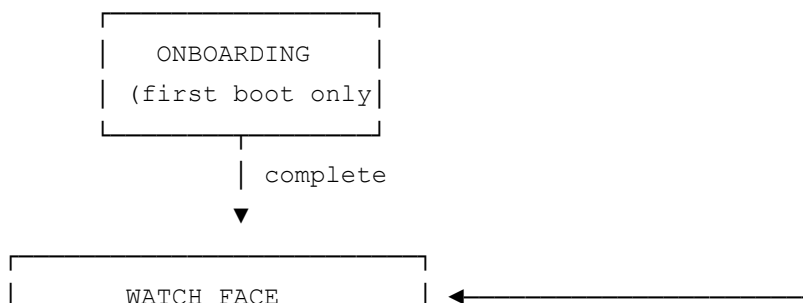

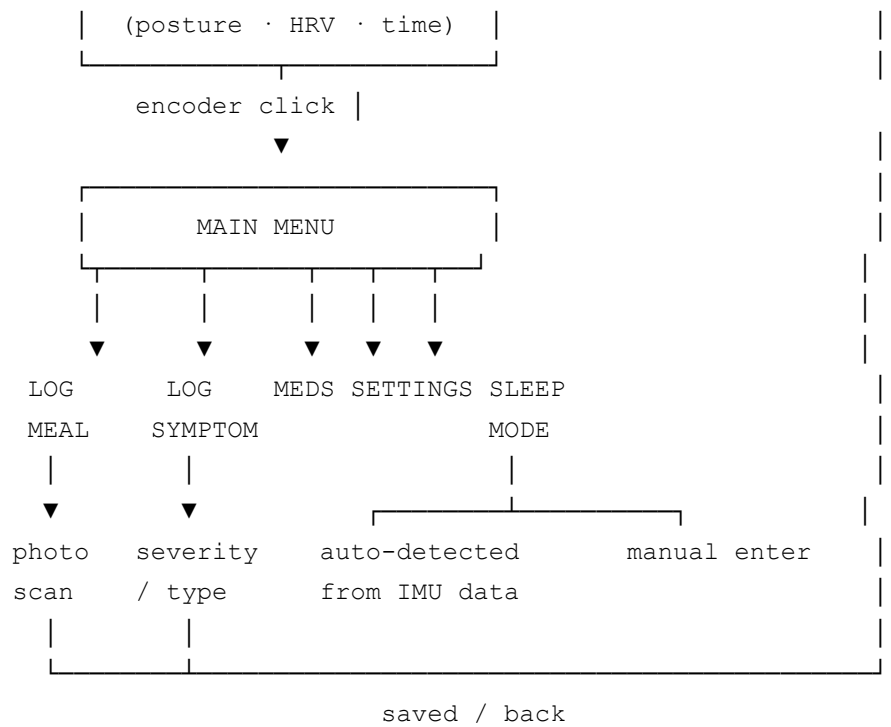

All OLED screens render at 128×64 px. Font: 1× and 2× scale system font. Invert highlight used for selected items. Screens auto-return to watch face after 30s inactivity.